

aiSon™ FOCUS

The world's only adaptive ultrasound pad.

✉ info@aisontechnologies.com

🌐 www.aisontechnologies.com

“ Well, it's the best ultrasound standoff pad that I have so far tested in my career and I have tested a lot and was never happy. ”

- Experienced sonographer & trainer of MSK sonography


● Better images with ease

The aiSon™ FOCUS' adaptive shape ensures continuous contact with the anatomy and a full field of view during the whole procedure. The pad moves effortlessly with the probe, always stays in place and increases image contrast for clearer visibility of structures.



● Reduced pressure

The aiSon™ FOCUS avoids the compression of the anatomy that could alter the diagnosis and reduces pain when scanning sensitive areas such as fractures or nerve lesions.

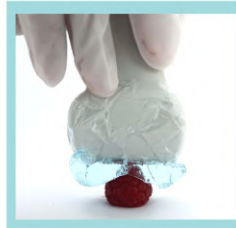

● Significantly less gel

The aiSon™ FOCUS increases patients' comfort while gaining improved quality images. The pad also lessens frequent re-gelling resulting in an uninterrupted workflow, and reduces the use of gel to bridge uneven surfaces or as a standoff.



● Re-usable & disinfectable

The aiSon™ FOCUS is usable up to 12 hours on multiple patients, protects the ultrasound sensor from disinfectant, and reduces maintenance costs.



● Easier interventions

The aiSon™ FOCUS provides significantly better view in semi-sterile guided interventions where disinfectant is used instead of gel. When sterile conditions are required, a sterile probe cover can be placed on top of the aiSon™ FOCUS easily.



● Versatile applications

The aiSon™ FOCUS works with all customary ultrasound probes and can be used either as an adaptive pad or as a standoff pad. The pad is beneficial for innumerable specialties such as radiology, gynecology, angiology, anesthesiology, pain medicine, pediatrics, rheumatology, dermatology, emergency medicine, hand & plastic surgery, ENT, ophthalmology, veterinary medicine, etc.



Product Information

Type	Non-sterile, consumable, reusable for 12hrs
Regulatory status	Class 1 (EU MDR) CE conform
Compatibility	Works with all customary ultrasound probes
Shelf life	24 months
Storage	Room temperature
Quantity & packaging	50 units/box
Weight	10g/unit
Disinfection	Disinfectant without pyridine
Ingredients	Recyclable plastic & water-based filling, free of latex & phthalates, disposal into customary recycling bin
GTIN code	7649988699412
Owner & manufacturer	Aison™ Technologies AG Wiesenstrasse 10A 8952 Schlieren Switzerland
EU authorized representative	Aison™ Technologies Poland Sp.z o.o. Gęsia 8 31-535 Krakow Poland

